


Image Manager or Image & Video Manager

Digital Storage

Fast – Up to date – Everywhere

Our vision of the connected healthcare is that medical professionals should have access to high quality image and video related diagnostic services of their choice, at each moment and each place in a secure, intuitive and easy method.

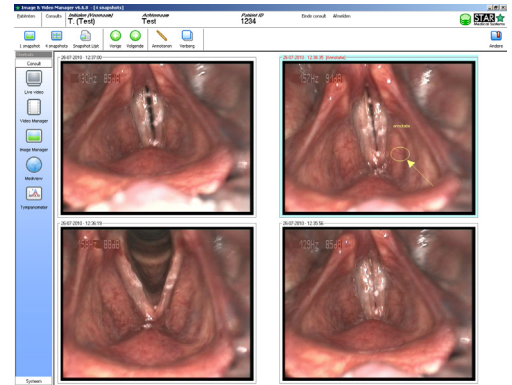

The complete image archive within reach:

- The Image Manager is a software application that allows images to be recorded and archived.
- The Image & Video Manager is a software application that allows video and images to be recorded and archived.
- Patient information can be easily linked with the images and video clips recorder with the Image (and Video) manager made .
- The Image (and Video) Manager can be used for all specialties in which 'visual-light' should be included.

Linking possibilities:

STAR Medical Systems offers hospitals and physicians a complete solution to enable links with the Hospital Information System (HIS), Picture Archiving and Communication System (PACS) and Electronic Patient Record (EPR).

- Video and images can be linked to the EPR through the STAR HL7 Gateway or STAR DICOM Gateway.
- The name and address of the patient from the hospital information system can be acquired on the Image (and Video) Manager workstation using the STAR HL7 link.
- Worklist are with DICOM worklist link easy to use on the Image (& Video) Manager workstation.
- STAR Medical Systems offers a complete solution for storage of 'visual-light' images in PACS (II).
- Acquiring images from different locations is possible with the WebIVM. The images can be accessed anywhere and anytime. This could be your own workstation, an operating room or any other place with a network connection.
- By linking to our Hygiene Tracker software, the status of the different scopes are recorded by department- or hospital level.
- Using our Clinical Reporter structured reporting can be done , enclosing relevant images from the Image (and Video) Manager.
- With one touch you can, at patient level, export all (video) images to CD / DVD, USB stick or other external data carriers.

Optimization ease of use

This is just a brief summary of the possibilities that the Image (and Video) Manager provides. The Image (and Video) Manager is especially designed for ease of use, flexibility and wide user options

Cut and Paste.


The Image & Video Manager allows you to capture and store images and video streams. Video streams can be viewed frame by frame and edited. You can also cut single frames.

Annotation

You have the possibility to make annotations and the possibility to make the annotations visible or not on the original image. Edited images can be exported for use in presentations or educational material.

Accurate overview of image (and video) reports

The Image & Video Manager shows multiple images on a single screen. This allows you to compare images from the patient history with the actual situation on the same screen.

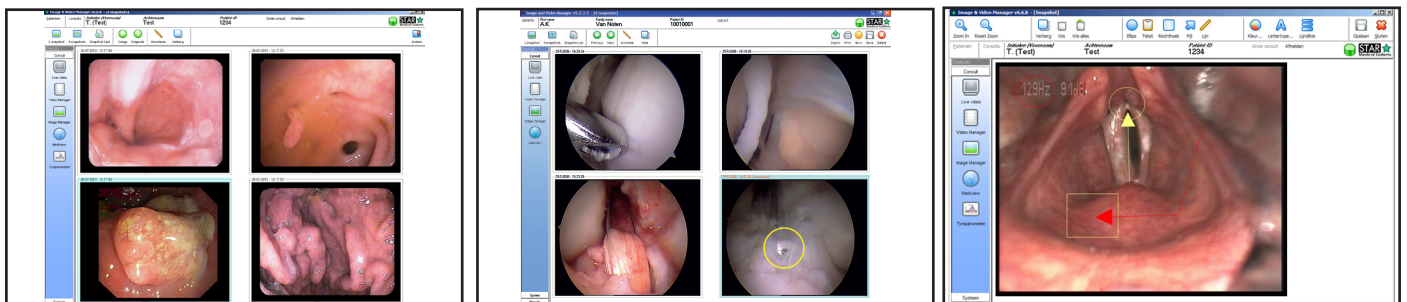

Especially for specialist

departments; Gastroenterology

- The Image (and Video) Manager can be controlled using the buttons on your endoscope.
- There is a link possible to our digital reporting software (Clinical Reporter) so the images stored with the Image (and Video) manager easily can be integrated

Orthopaedics

- With a headset, audio commentary can automatically be added to the video recording.

Ear Nose and Throat

- Capturing images and / or video using the foot pedal
- Full integration with the Entermed Futurent treatment units
- Link with relevant ENT measuring equipment